

WelComb® has all you need and nothing you don't!

NO FEAR. NO DRAMA. NO PESTICIDES. NO GAPS.


To get rid of lice and nits, you need the right comb. Most of the combs on the market slip right past the .3mm-wide nits (eggs). With rigid teeth just the right distance apart and a beveled edge that gets right down to the scalp, the WelComb® always gets the job done!

The comb nurses can confidently recommend!

Other Combs vs. WelComb®

SEE THOSE GAPS?

They can leave lice and nits behind.

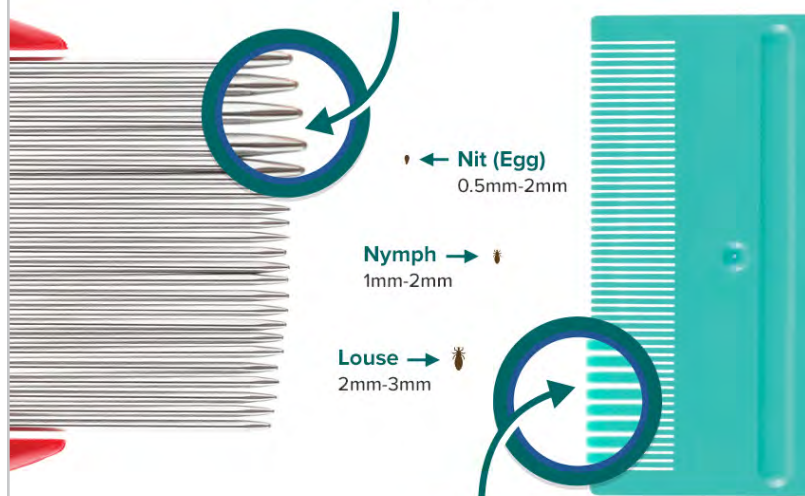

With WelComb®,
NO GAPS!

WELCOMB®
LICE AND NIT REMOVAL COMB



SHOP FOR WELCOMB® AT

DOLLAR GENERAL®

OR VISIT

www.welcomb.com

What Nurses Need to Know About the Wet Combing Method

Did you know? It's the **comb**, not the chemicals, that gets rid of lice and nits. Wet combing is all about removing the live lice and the viable nits from the head. Simply removing the live lice is enough to kill them and removing the viable nits prevents them from hatching.

Remember, there is **no** pesticide or natural solution that can harm or kill nits. That means the **only** way to ensure that nits won't hatch and continue to spread is to physically remove them. And while it is low cost and safe, the wet combing method takes time and patience.

For the wet combing method, all you need is:



Full wet combing instructions, including videos in English and Spanish, can be found at welcomb.com/instructions and welcombspanol.com/instrucciones



Water because lice freeze when wet, making it easier to capture and remove them.



Conditioner to lubricate the hair and make it much easier to comb through.



WelComb® to physically comb out lice and nits. With a gap between the teeth that's .15mm, it's small enough to capture even the smallest nit.

SHOP FOR WELCOMB® AT

www.welcomb.com

MacGILL®
SCHOOL NURSE SUPPLIES

WELCOMB®
LICE AND NIT REMOVAL COMB

