


Prepared For:
Jinyu Li
Tianjin TEDA Co. Ltd.
South Yinhe Rd., Beichen Dist.
Tianjin City 300400
P.R. CHINA

Submitted By:
Nelson Laboratories, Inc.
6280 South Redwood Road
Salt Lake City, UT 84123-6600
801-290-7500
Page 1 of 2


**BACTERIAL FILTRATION EFFICIENCY
AND DIFFERENTIAL PRESSURE – FINAL REPORT**

Laboratory Number:	359507
Procedure Number:	SOP/ARO/007L.1
Sample Source:	Tianjin TEDA Co. Ltd.
Sample Identification:	Refer to Table 1
Deviations:	None
Statement of Uncertainty:	If applicable, available upon request
Andersen Sampler Flow Rate:	28.3 L/min. (1 CFM)
BFE Conditioning:	4 hours minimum at 21 ± 5°C and 85 ± 5% relative humidity
Sample Received Date:	15 Jan 2007
Lab Phase Start Date:	16 Jan 2007
Lab Phase Completion Date:	28 Mar 2007
Report Issue Date:	29 Mar 2007
Results:	Refer to Table 1

The Bacterial Filtration Efficiency (BFE) procedure is performed to determine the filtration efficiency of various materials and filtration devices using a challenge organism of *Staphylococcus aureus*. This procedure complies with ASTM F2101. The Differential Pressure (Delta P or ΔP) test is performed to determine the air exchange differential (breathability) of porous materials.



Technical Reviewer

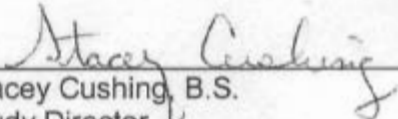

Stacey Cushing, B.S.
Study Director

29 Mar 2007

Study Completion Date

bav



TABLE 1. Results

SAMPLE IDENTIFICATION	ΔP (mm H ₂ O/cm ²)	PERCENT BFE
TDK-25 - 1	2.0	>99.9%*
TDK-25 - 2	1.7	99.9%
TDK-25 - 3	1.8	>99.9%*
TDK-25 - 4	2.1	>99.9%*
TDK-25 - 5	2.2	>99.9%*

CONTROL AVERAGE: 2235 CFU

MEAN PARTICLE SIZE (MPS): 3.2 μm

* There were no detected colonies on any of the Andersen sampler plates for this sample.